

Comfort Duett®

The compact multipurpose personal amplifier



Bringing you closer together

The compact multipurpose personal amplifier.

Comfort Duett® is a small and easy-to-use, rechargeable assistive listening device that provides superior sound. It is stylish and ergonomically designed. The Comfort Duett® amplifies speech and sounds from a variety of sources like TV, radio, telephone, etc.

The innovative base unit serves as an overnight charger and as an interface to audio equipment such as your TV and telephone. When you use earphones, the Comfort Duett®, is very discreet and resembles a stylish music player.

Comfort Duett can be used with headphones, earbuds or a neckloop for those who use hearing aids with a built in T-coil.

All this makes Comfort Duett® the ideal communication tool for a majority of users.



Features

-  In small groups
-  Watching TV
-  Talking on the phone



Comfort Audio
HEAR THE FUTURE

In everyday situations




Personal Amplifier



Talking on the phone



At a meeting

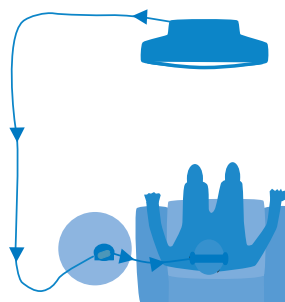

Watching TV

Technical facts

Operating Time:	Up to 20 hours, depending on volume level
Charge time:	1 hour charge = 1 hour use; (14 hours from empty to fully charged)
Weight:	1.75 oz., including battery
Battery:	Rechargeable NiMH
Size:	3.38 x 1.50 x 1.13 inches
Operating spectrum:	275 - 9000 Hz.
Max. amplification:	60 dB, with max gain and bass reduction controls
Max. SPL:	133 dB
Color:	Metallic blue
Warranty:	2 years limited manufacturer's warranty
Built-in T-coil:	Can be used in public areas that have induction loop system installed